‘Let every man judge according to his own standards,
by what he has himself read, not by what others tell him.’
Albert Einstein (1949) *The World As I See It*

ZSI
Zephyr Surgical Implants

SWISS MADE

- ***ZSI** is the **only European manufacturer** of artificial urinary sphincters and of malleable and inflatable penile implants, as solutions for male urological problems (incontinence), sexual complications (Erectile Dysfunction), and gender identity procedures (prostheses for trans men and women).*

MADE IN SWITZERLAND
In the heart of Europe
ZSI, Zephyr Surgical Implants
Route des jeunes 4bis,
Les acacias-Geneva
CH-1227, Switzerland

PRODUCTS

The hydraulic / inflatable penile implant in three components (PI) ZSI 475 is a solution for problems related to erectile dysfunction and impotence.

A malleable / semirigid penile implant for men. ZSI 100 is a solution for problems related to erectile dysfunction and impotence.

ZSI offers solutions adapted for plastic surgery to satisfy the specific demands of transsexual / transgender men and women.

A hydraulic / inflatable penile implant of three components for trans men ZSI 475 FTM

A malleable / semirigid penile implant for trans men ZSI 100 FTM


A malleable / semirigid implant for metoidioplasty for trans men ZSI 100 D4

ZSI
Zephyr Surgical Implants

50% of men over 40 have some degree of ED*.

Erectile dysfunction causes:

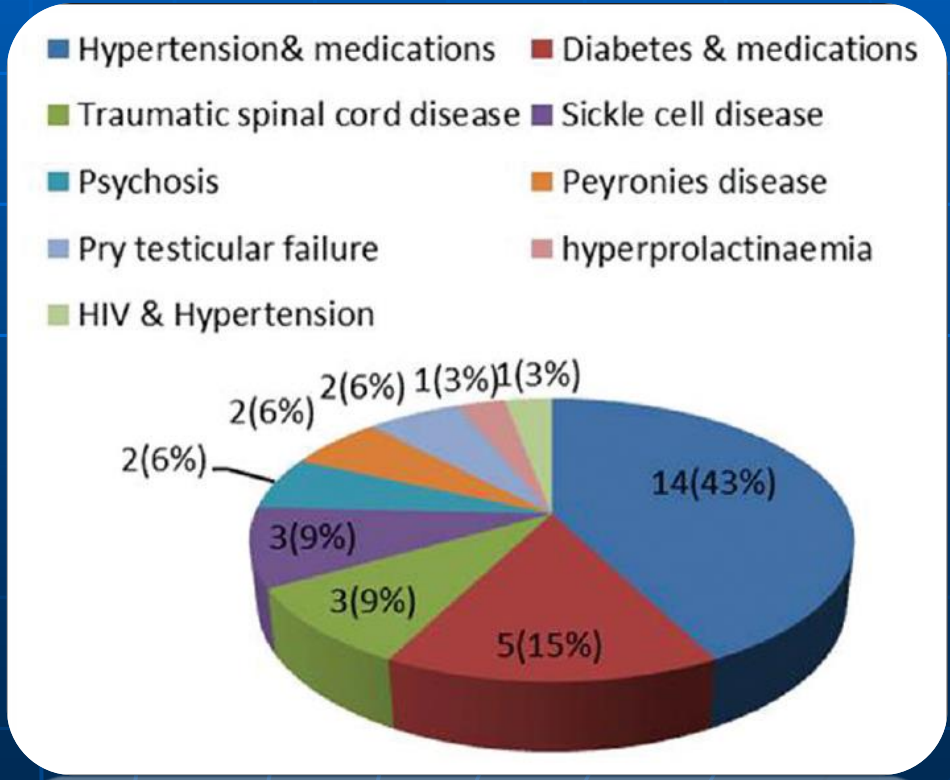ED treatment strategy

Always the same what ever the cause:

- Medical history for psychogenic or organic etiology
- Physical exam looking for blood pressure in cavernous artery or corpus cavernosa fibrosis
- Biological tests: Blood glucose, LH-FSH, Prolactin, Testosterone (not useful main of time)

TREATMENT

- Psychogenic: Oral therapy/per-os
- Organic (diabetes, radical prostatectomy, colectomy, vascular surgery….)

always the same strategy:

- ✓ Oral therapy/ per os or Vacuum
 In case of failure
- ✓ Intra cavernous injection
 In case failure
- ✓ Penile implant

ZSI
Zephyr Surgical Implants

ORAL THERAPY/PER-OS:

- VIAGRA (Sinldenafil)
- CIALIS (Tadalafil)
- LEVITRA (Vardenafil)
- SPEDRA (Avanafil)


VACUUM


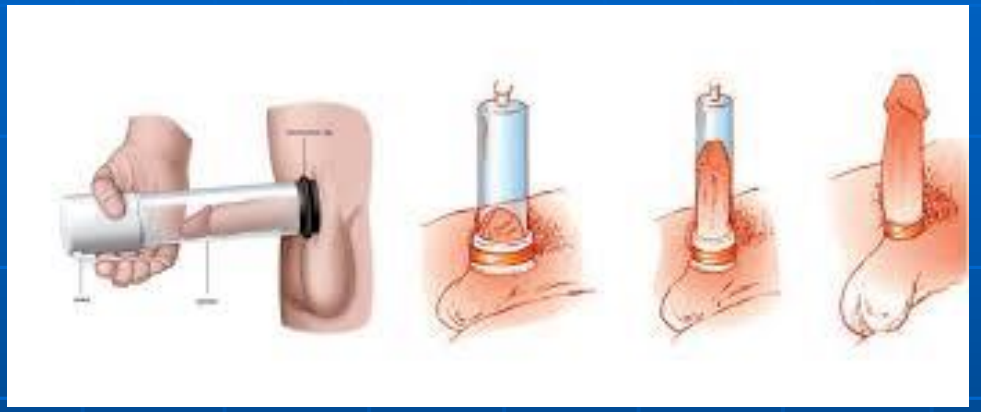
INTRA-CAVERNOUS INJECTION

Corpus cavernosum


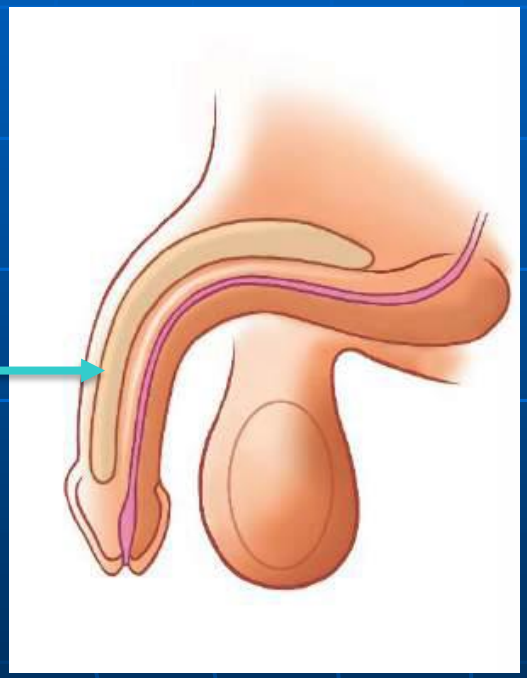
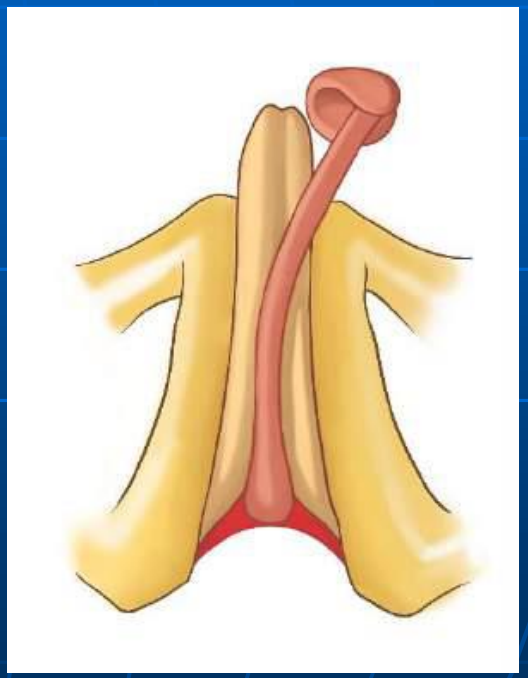
ZSI
Zephyr Surgical Implants

PENILE IMPLANT

- They are intra-cavernous implant.
- Glans is respected, it keep its sensibility.
- If the patient still have his prostate he can ejaculate

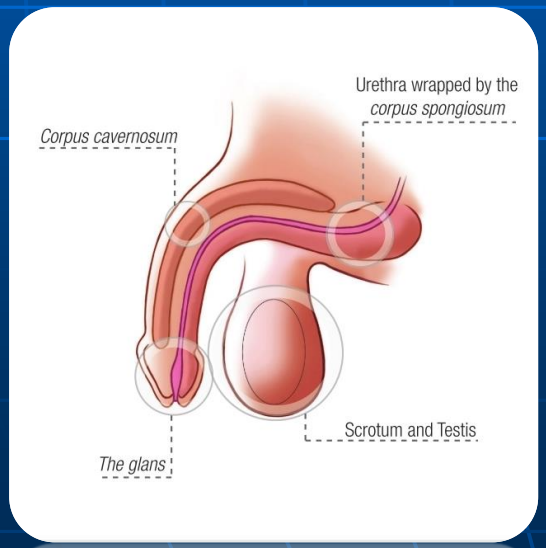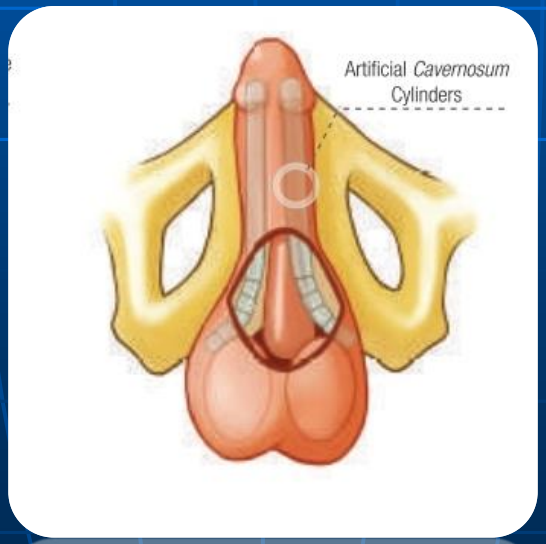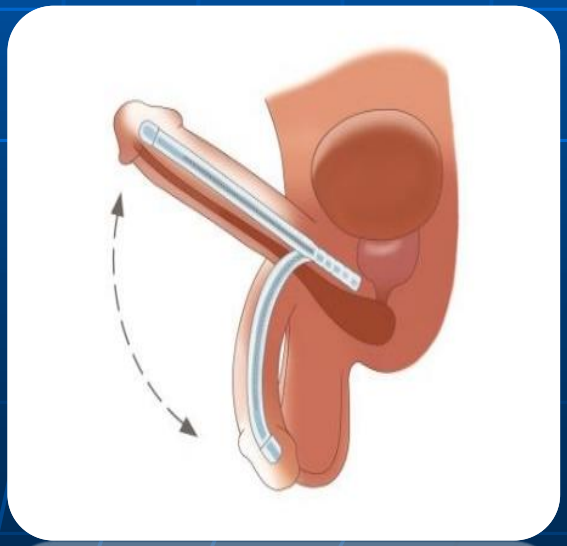SEMI-RIGID / MALLEABLE PENILE IMPLANT

ZSI
Zephyr Surgical Implants


- The cylinders are inserted within the corpus cavernosum, residing on either side of the penis
- No tissue is removed or created to insert the cylinders
- They are placed within an already existing space of the penis, where previously blood would flood and cause the penis to be erect

The ZSI 100 does not affect the:


- Patients capacity to urinate or ejaculate
- Natural sensitivity of the penis nor of the Glans


SEMI-RIGID / MALLEABLE PENILE IMPLANT FUNCTIONING

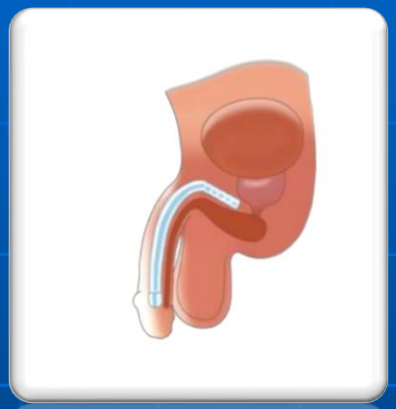
flaccid penis

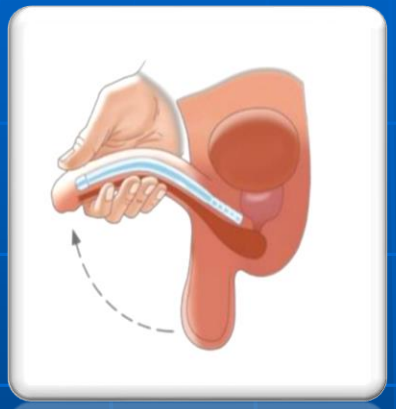With hand straighten the penis from bent position

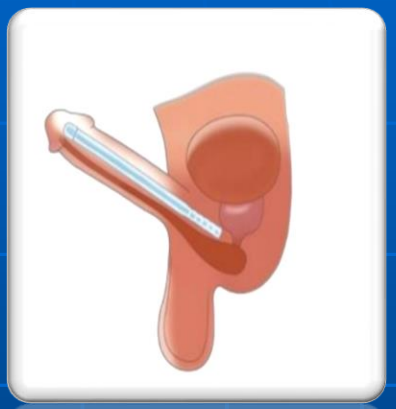Erect penis

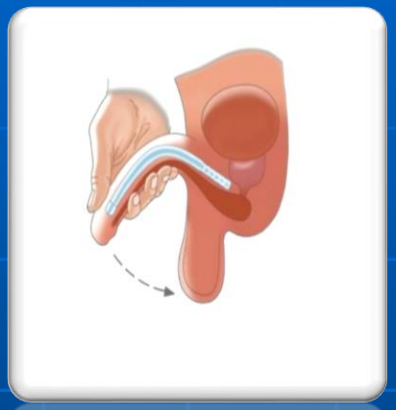With a hand return the penis to bent position

The device is implanted in the corpus cavernosum.

ZSI
Zephyr Surgical Implants

MALLEABLE PENILE IMPLANT

ZSI 100


A Solution for erectile dysfunction

MALLEABLE PENILE IMPLANT
ZSI 100

ZSI
Zephyr Surgical Implants


Manufactured at an optimum density for both durability and malleability, providing a lifelike erection.

Reinforced Nitinol (Nickel-Titanium) rod interior to provide best long term use.

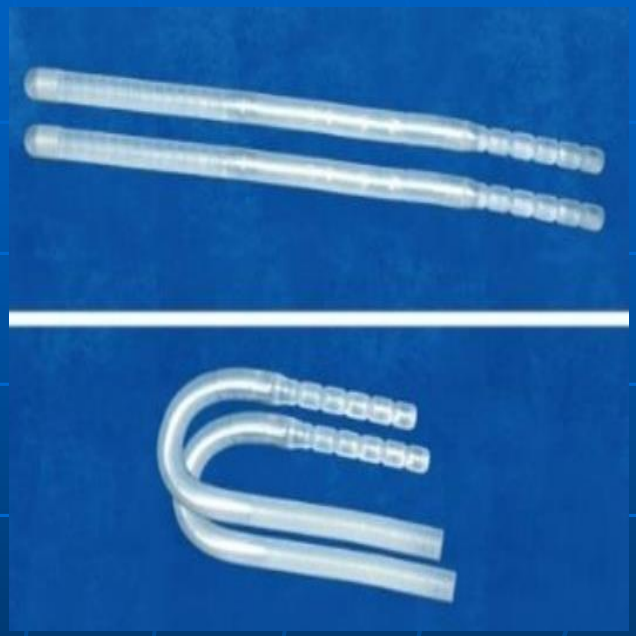ZSI 100


- Length is 25 cm. It can be reduced up to 12 cm
- It can be adjusted from both distal and proximal part
- Always cut distal part first and then proximal part if necessary. This is important for cable position in corpus cavernosum.
- 3 diameter are available: 9mm, 11mm, 13mm
- In case of mistake, extensors are provided.

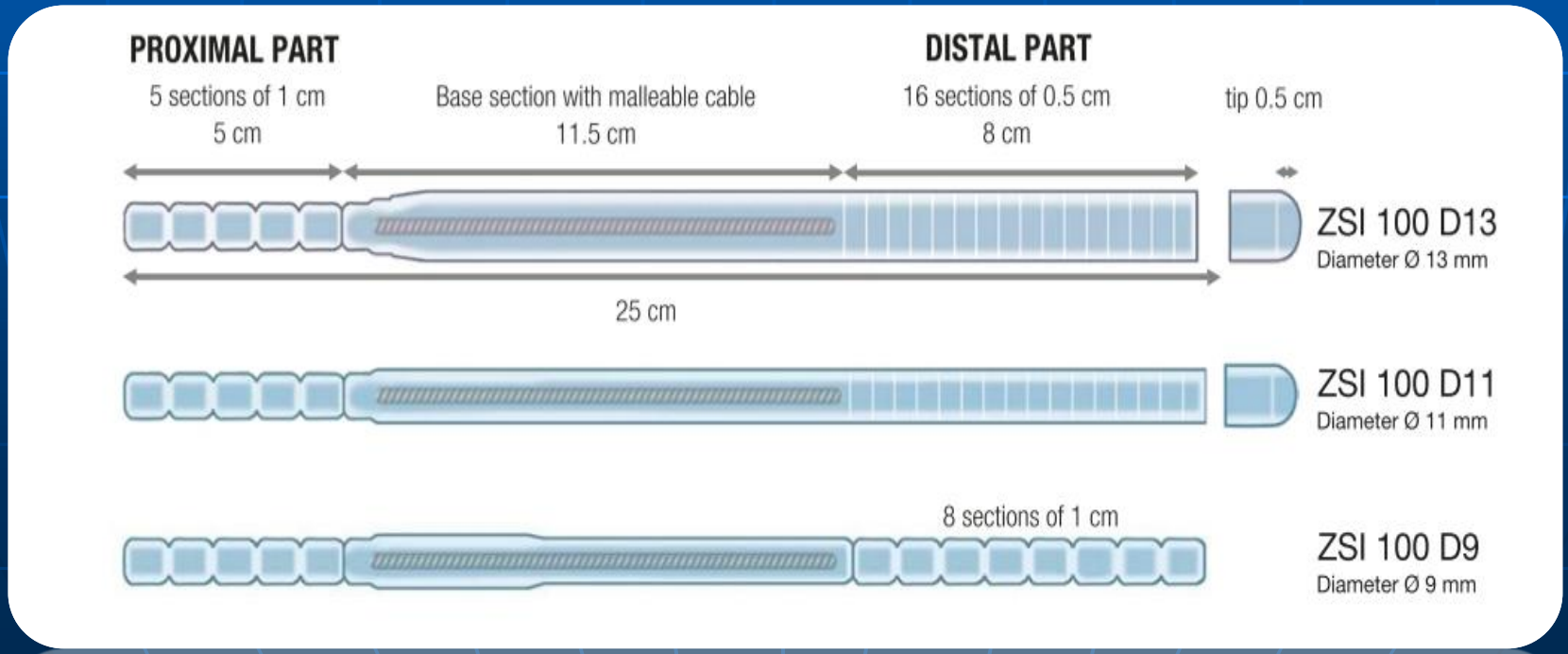
ZSI 100

ADJUCTING FROM BOTH DISTAL & PROXIMAL PART

The surgeon can reduce the length from the distal or/and proximal part, to suit the individual patient's anatomy or expectations. This allows the ZSI 100 device to be modified to maximize patient comfort.

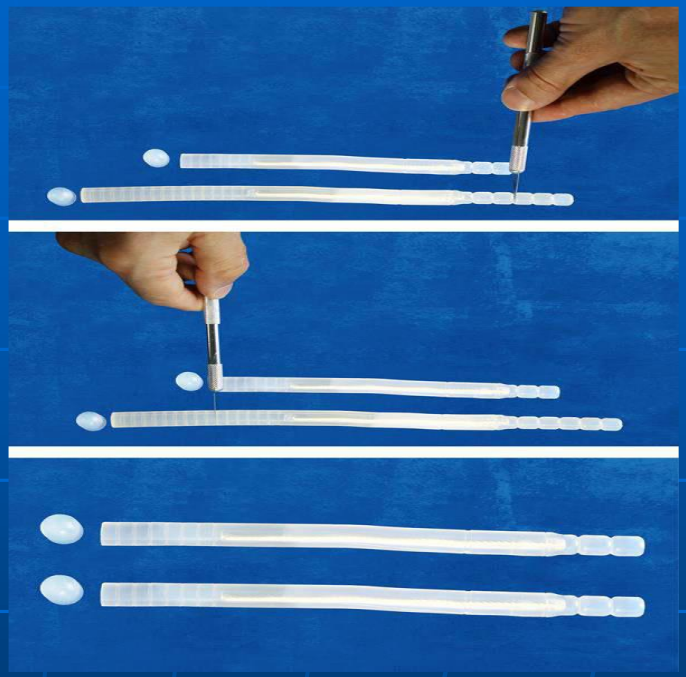
- Extensors are available in the kit only to repair an excess cut of the cylinders.

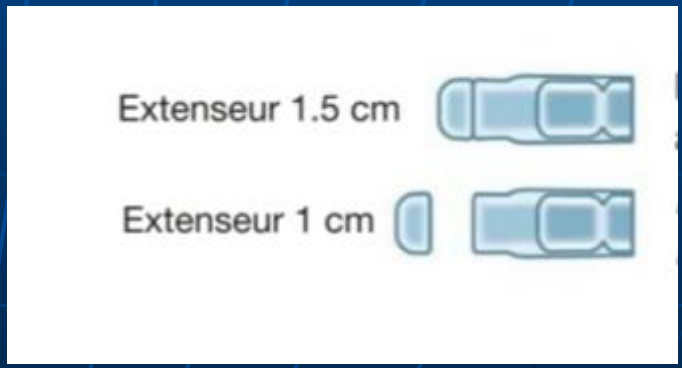ZSI 100

ZSI
Zephyr Surgical Implants

PROTECTION OF THE ENTIRE PENIS

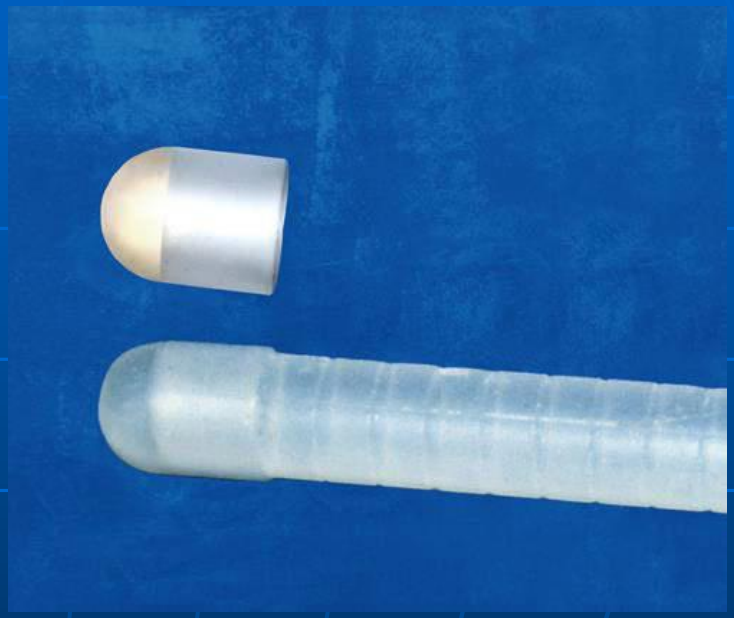
The ergonomically shaped stopper provides a secure connection to the gland after cutting at the right length.
NUSIL is the most secure biocompatible silicone that reduces the risk of infection and erosion.

ZSI 100 Advantages

- **Excellent memory form**

1) The ZSI 100's excellent shape memory is an important feature which increases patient comfort. The less readily the implant resumes its shape after bending, the more comfortable it will feel in a resting position, in the patient's trousers.

2) Implants without good shape memory will be uncomfortable for everyday use.

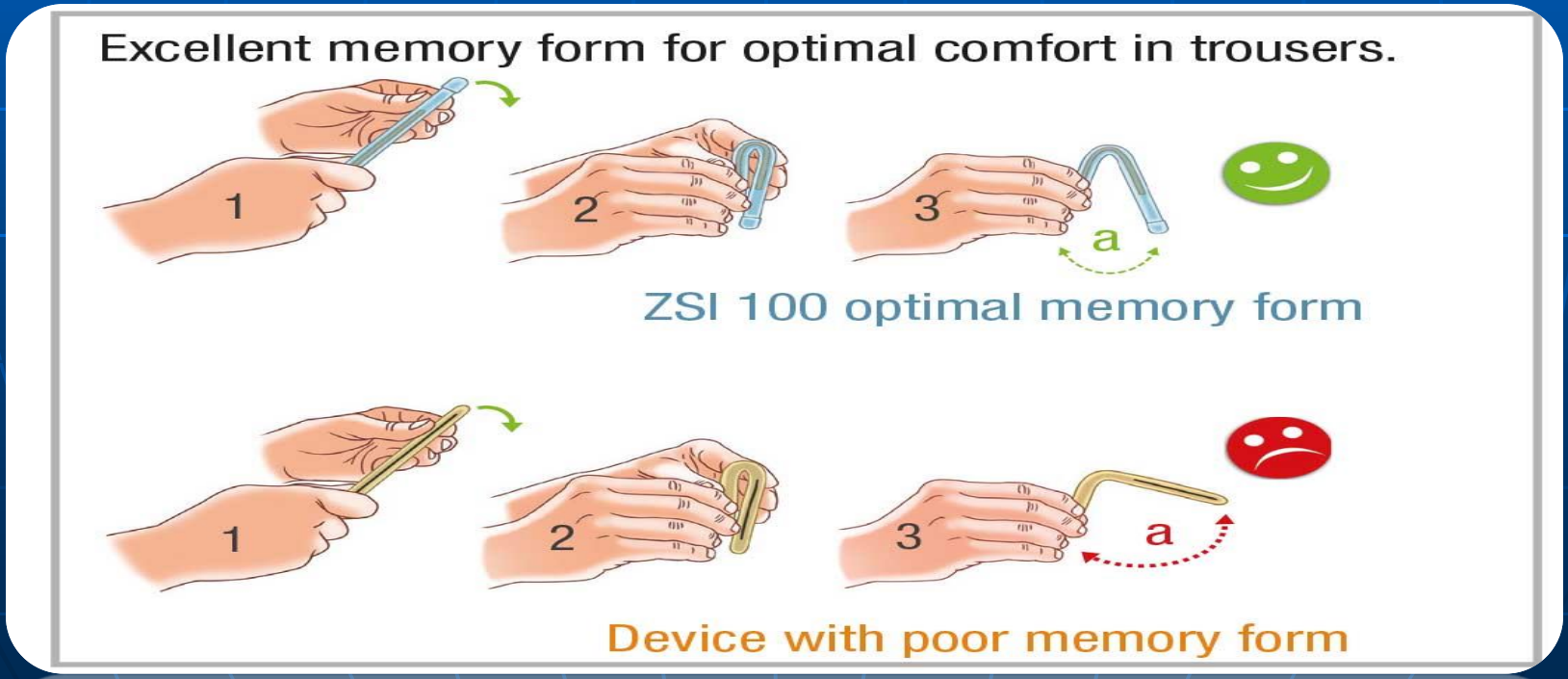
ZSI 100 Advantages

- **FLEXIBLE DISTAL PART**

1) Pressing the side of a penis implanted with the ZSI 100 results in a natural bending motion which feels like that of a natural penis.

➢ **LACK OF FLEXIBLE DISTAL PART**

2) Pressing the side of a penis implanted with an overly rigid distal part does not result in a natural bending motion. It will feel unnatural.

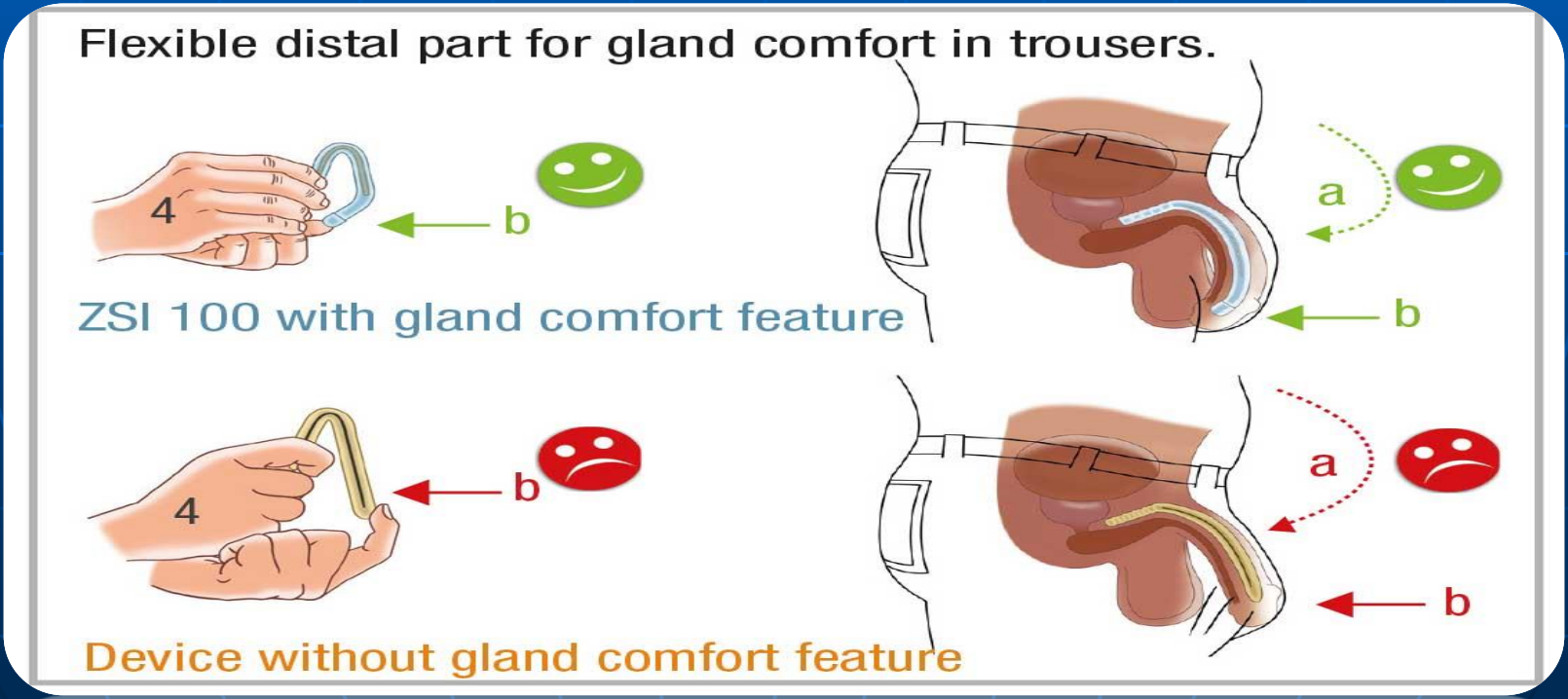
ZSI 100 Advantages

- **RIGIDITY FOR INTERCOURSE**

1) 2 cylinders of ZSI 100 provide the ideal rigidity for intercourse. Any size ZSI 100, from 12 to 25cm, will provide satisfactory results for intercourse.

➢ **LACK OF RIGIDITY FOR INTERCOURSE**

2) device lacking rigidity might be very comfortable in a resting position, but would be useless for intercourse.

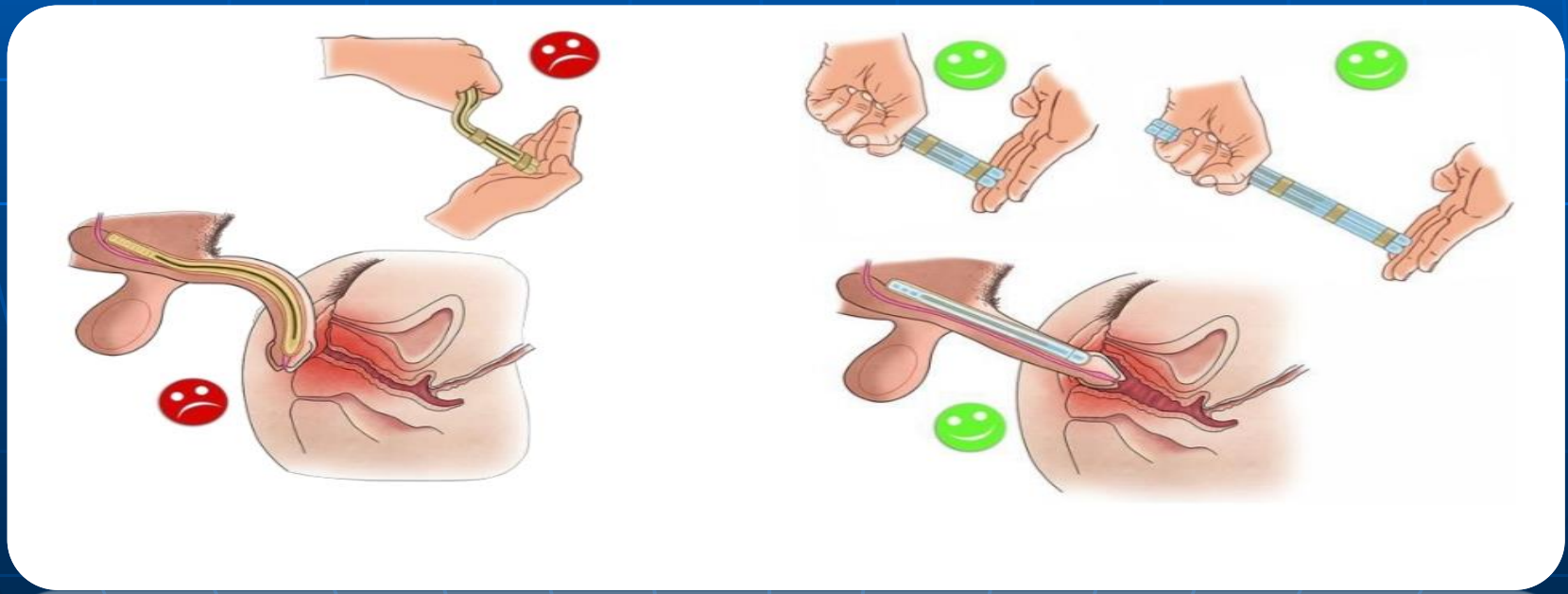
ZSI 100 Advantages

- **GLANS COMFORT**

1) The ZSI 100 has a flexible distal part, allowing for optimum gland comfort in a flaccid position


➢ **LACK OF GLAND COMFORT**

2) A lack of flexibility in the distal part of the device will produce excess rigidity, leading to gland discomfort in the resting position

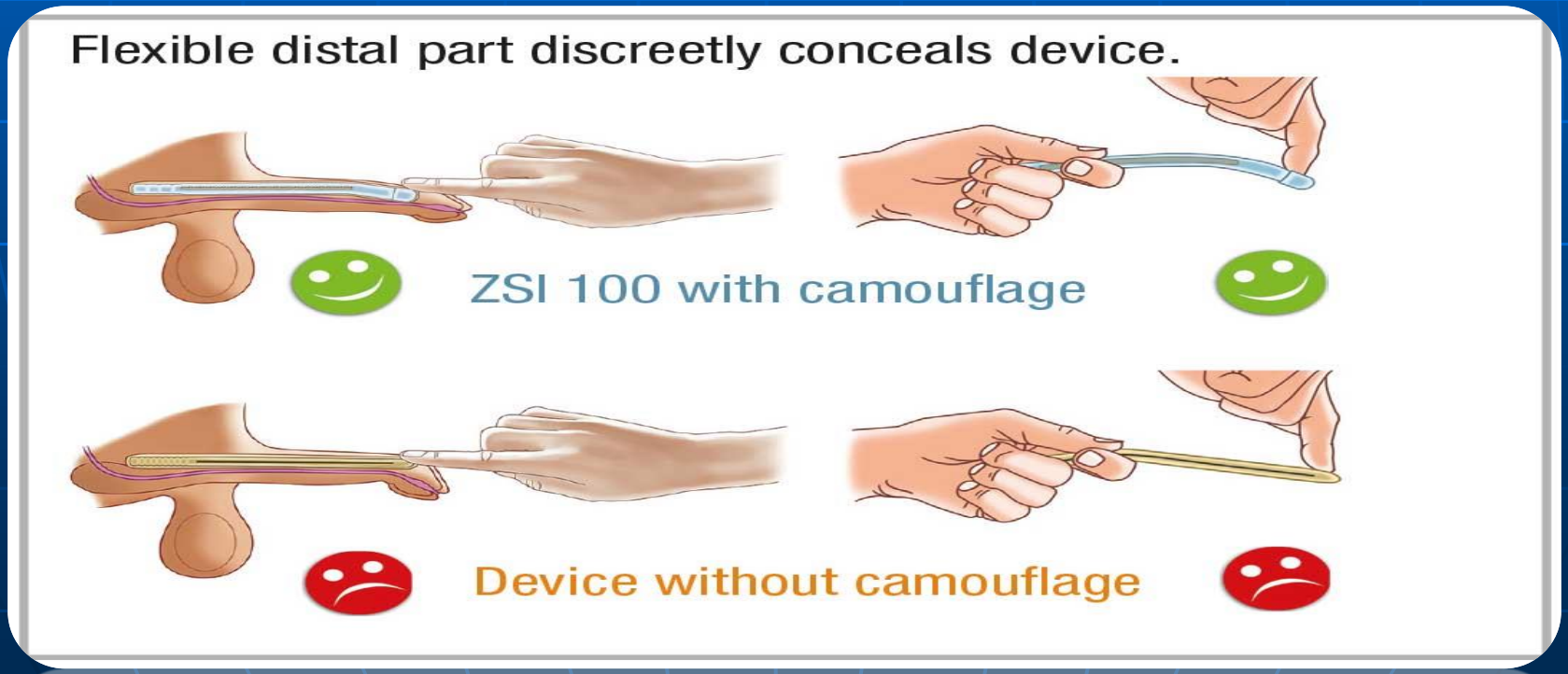
ZSI 100 Advantages

- ✓ Pair of implants, fits all lengths,
- ✓ Adjustable proximal and distal from 12 to 25 cm,
- ✓ Easy implantation,
- ✓ Optimum patient comfort,
- ✓ Excellent memory form,
- ✓ Good endurance,
- ✓ Good rigidity,
- ✓ Good plasticity,
- ✓ Comfortable in the rest position,
- ✓ Effective in erection.
- ✓ Safe in MRI & airport scanners

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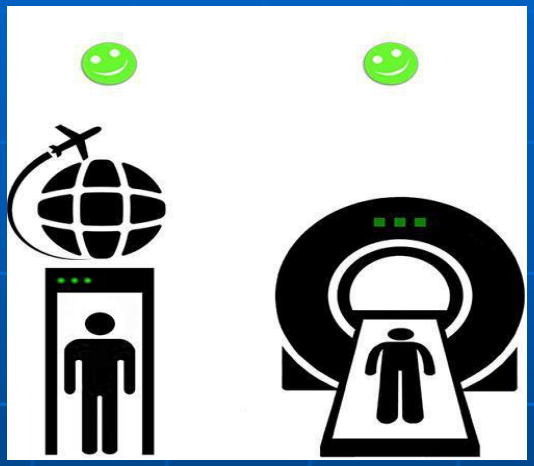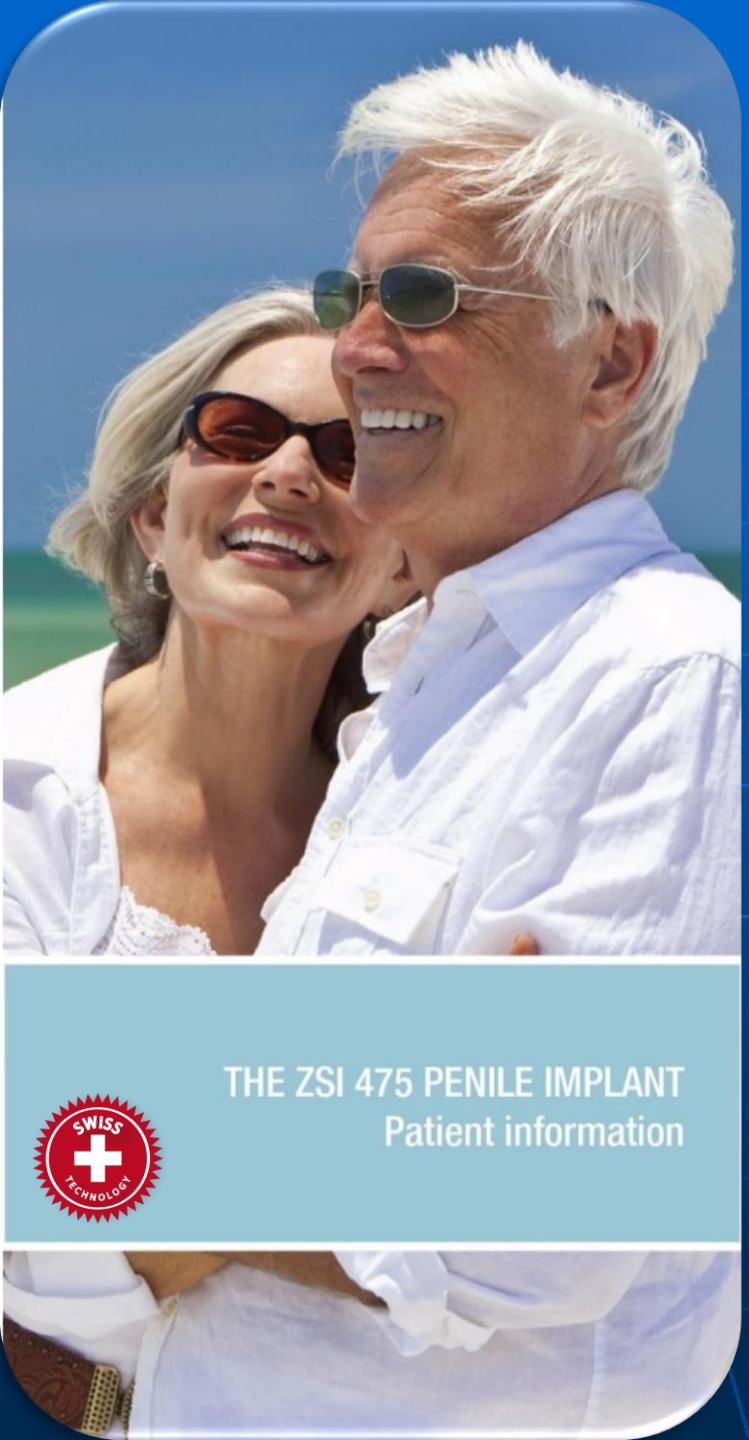INFORMEDCONSENT, main information

- ✓ Erection of corpus cavernous but no glans erection
- ✓ Glans keep its sensitivity
- ✓ Penis size the same before and after the procedure
- ✓ <10% infection, <10% erosion but there is a risk the patient must accept before the procedure.
- ✓ Risk of perineal pain for few weeks.

ZSI
Zephyr Surgical Implants

THANKS FOR YOUR ATTENTION

THE END