
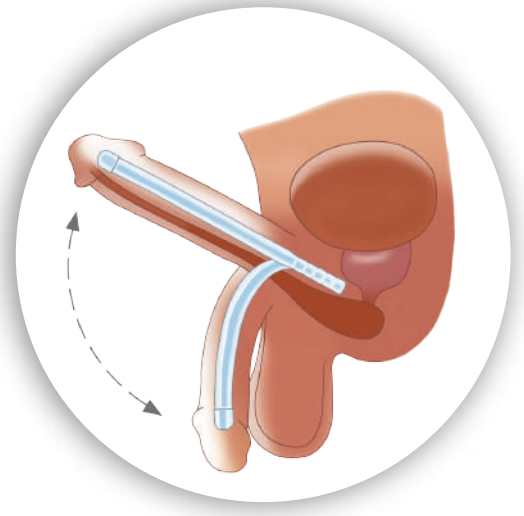


ZSI 100

MALLEABLE PENILE IMPLANT

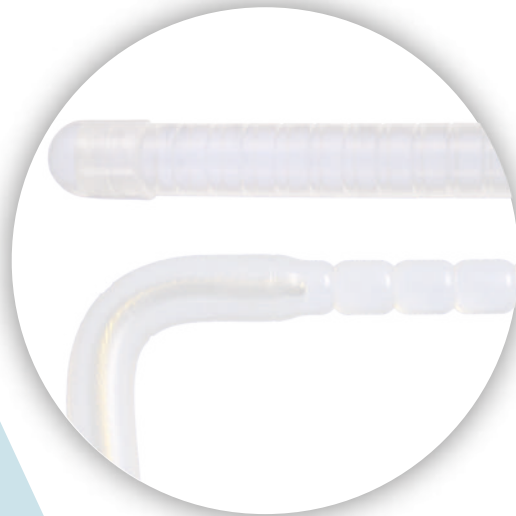

A solution to erectile dysfunction

- pair of implants, one size fits all.
- for more information, visit www.zsimplants.ch

ZSI is ISO 13485 : 2004 / ZSI 100 is EC certificated

Local Distributor

ZSI (Zephyr Surgical Implant Sàrl)
C/Ulmann, route des Jeunes 4bis
1227 Les Acacias - Switzerland


Flyer, English version
Version 4, September 2016



THE ZSI 100 SEMI-RIGID PENILE IMPLANT DESIGNED TO PROVIDE EXCELLENT COSMETIC RESULTS

Sep. 2016 V4

ZSI 100 ADDED VALUE CAN BE TESTED:

Ask your local distributor to test our ZSI 100.

Excellent memory form for optimal comfort in trousers.

ZSI 100 optimal memory form

Device with poor memory form

Flexible distal part for gland comfort in trousers.

ZSI 100 with gland comfort feature

Device without gland comfort feature

Flexible distal part discreetly conceals device.

ZSI 100 with camouflage

Device without camouflage

Together, the cylinders can withstand up to 1kg pressure. Ideal rigidity for intercourse.

PAIR OF IMPLANTS, 11MM IN DIAMETER FITS ALL PENILE SIZES ADJUSTABLE PROXIMALLY AND DISTALLY FROM 12 TO 25CM FOR EASE OF IMPLANTATION AND OPTIMAL PATIENT COMFORT

DISTAL PART
5 mm by 5 mm.
Total length 80mm.
Adjustable by cutting.

CENTRAL PART
Total length 120 mm.
Bendable rod interior.

PROXIMAL PART
10 mm by 10 mm.
Total length 50mm.
Adjustable by cutting.



- 11 mm diameter fits all standard patients.
- Manufactured from NUSIL silicone for optimum durability.

The patient can manipulate the reinforced interior cable to produce an erection.

The ergonomically shaped stopper provides a secure connection to the gland.

